

Pelican Plus® Pool & Spa Hoist with Platform Lift



Designed and manufactured in Australia by Para Mobility, and tested in accordance with ASA 3581, the Pelican Plus Pool & Spa Hoist with Platform Lift offers direct access to pools via wheelchair without any manual lifting.

At Para Mobility, we understand that all pools are different, so we work with you to ensure the Pelican Plus Pool and Spa Hoist meets your specific pool needs. The unit can be mobile or fixed permanently at the edge of the pool or spa and is easily installed on the pool deck by a handyperson, or alternatively we can take care of it for you via our Pelican® Installation & Service.

FEATURES

- Complies with 'Disability (Access to Premises – buildings) Standards 2010 / D3.10 Swimming Pools'
- Easily installed (Pelican® Installation & Service is available)
- Easily wheeled from one area to another for multiple use and or storage (mobile unit)
- Platform attachment can be independently detached and stored
- Removable platform, interchangeable with a sling or aquatic wheelchair with the Pelican release adapter
- Fully user operated via hand control
- Standard lifting capacity 200kg
- Cover for platform attachment available (additional cost)

SPECIFICATIONS

- 1500 wide x 1100 deep
- Handrail height 900mm
- Post length dependent on pool depth up to 1.3m deep
- Safe Working Load: 200kg
- Batteries: Linak – 24v re-chargeable battery
- Power rotation 360 degrees
- Total Weight: Approx 90kg
- Floor Mounting Plate: 450 diameter
- Construction: Stainless steel, powder coated to colour of choice
- Hand control: Linak – push button
- Install: Anchored to pool deck with stainless steel anchors and eye bolts

OPTIONAL EXTRAS

- Interchangeable hanging bar with floating sling and wheelchair attachment

Our Pelican Plus Pool & Spa Access Hoist with Platform Lift ensures that pools with a perimeter of 70 meters or more comply with the Disability (Access to Premises – Buildings) Standards 2010. According to D3.10 Swimming Pools (c), pools with a perimeter over 70 meters must provide at least one accessible water entry/exit, such as a platform lift and an aquatic wheelchair.




Para Mobility
Making everyday living easier

Phone 1300 444 600

Pelican® Installation & Service
Delivery, installation + servicing Australia wide
enquiries@paramobility.com.au
www.paramobility.com.au

LEARN MORE

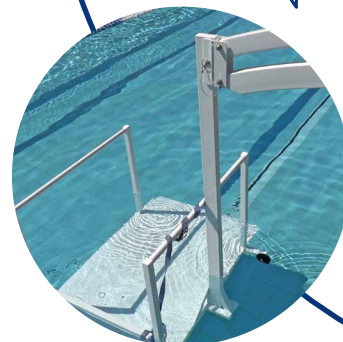

Freedom to explore